




SARS-CoV-2

Antigen Rapid Test Lollipop (Saliva)

A rapid test for the qualitative detection of SARS-CoV-2 nucleocapsid antigens in saliva specimens.

For professional *in vitro* diagnostic use only.

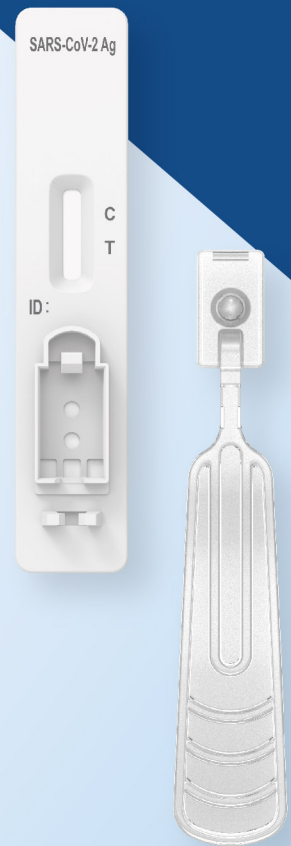

World Leader in Rapid Test Industry



Fast



Accurate



Easy to Use



Reliable



CE Marked



Flowflex SARS-CoV-2 Antigen Rapid Test Lollipop (Saliva)

The SARS-CoV-2 Antigen Rapid Test Lollipop (Saliva) is a lateral flow chromatographic immunoassay for the qualitative detection the nucleocapsid protein antigen from SARS-CoV-2 in saliva specimens directly from individuals who are suspected of COVID-19 by their healthcare provider within the first seven days of the onset of symptoms.

The SARS-CoV-2 Antigen Rapid Test Lollipop (Saliva) is intended for use by trained clinical laboratory personnel and individuals trained in point of care settings.

- **Specimen:** Saliva specimens
- **Test Time:** Results at 15 min.
- **Shelf life:** 12 months
- **Storage temperature:** 2-30°C



1T



5T



25T

Clinical Performance

Saliva Lollipop Specimens

The performance of SARS-CoV-2 Antigen Rapid Test Lollipop (Saliva) was established with 560 saliva Lollipop Specimens collected from individual patients. The results show that the relative sensitivity and the relative specificity are as follows:

Method	Results	RT-PCR (Nasopharyngeal Swab Specimens)		Total Results
		Negative	Positive	
SARS-CoV-2 Antigen Rapid Test Lollipop (Saliva) (saliva Lollipop Specimens)	Negative	450	20	470
	Positive	0	90	90
Total Results		450	110	560

Relative Sensitivity: 81.8% (73.3%-88.5%)*
Accuracy: 96.4% (94.5%-97.8%)*


Relative Specificity: 100.0% (99.2%-100.0%)*
*95% Confidence Intervals

Materials Provided


- Test Cassettes
- Saliva Lollipop
- Package Insert

Test Procedure

1. Test Steps

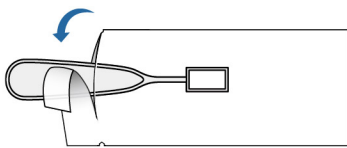


- 1 Open the saliva lollipop packaging at stick end.

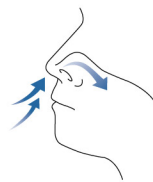


- 2 Drain mucus from the back of the nose and throat at least 3 times.



- 3 Take a deep breath in and cough hard, do not spit.

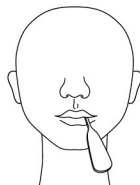


- 4 Insert a saliva lollipop into the mouth at least 30 seconds.

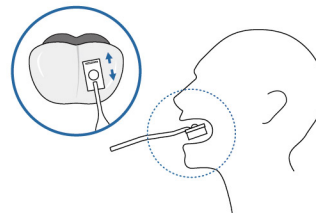

- 5 Wipe both **upper and down surface** of the tongue.

2. Running the Test

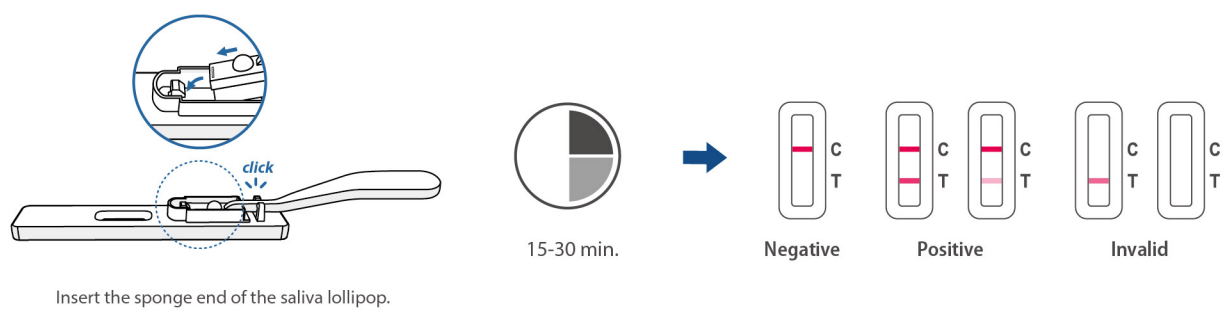

Ordering Information

Product Name	Catalog No.	Format	Specimen	Package
Flowflex SARS-CoV-2 Antigen Rapid Test Lollipop (Saliva)	L031-125R5 √	Cassette	saliva specimens	1 Test/Kit
	L031-125T5 √	Cassette	saliva specimens	5 Tests/Kit
	L031-125U5 √	Cassette	saliva specimens	25 Tests/Kit

√ CE Marked



ACON Biotech (Hangzhou) Co., Ltd.
No.210 Zhenzhong Road, West Lake District, Hangzhou, P.R.China, 310030
Tel: +86-571-8118 9706 / Email: information@aconlab.com.cn
www.aconbio.com